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Instruction manual

Manual de instrucciones

Mode d'emploi


EM EasyTwist™

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with bottle attachment

con adaptador para frascos

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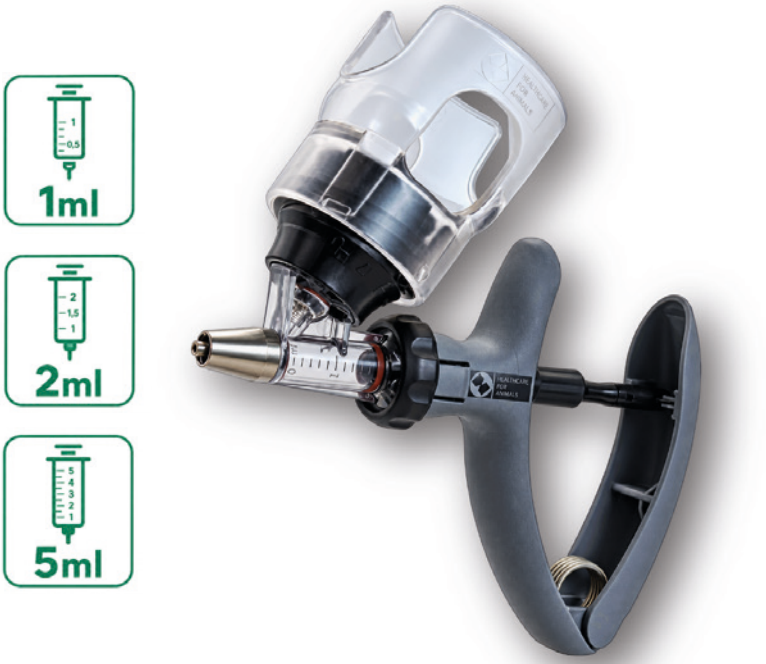
EM EasyTwist™
by HC4A

- Automatic self-filling syringe for injection
- with bottle attachment

A) ergonomically designed handle
B) plastic cylinder with bottle attachment
C) durable plastic easyMount basket
D) metal needle fitting, Luer Lock
E) dose adjustment sleeve
F) clip nut

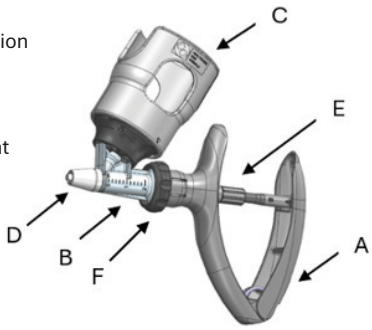INSTRUCTIONS FOR USE

The EM EasyTwist™ with bottle attachment is designed for daily use and longterm durability. Wear parts such as the O-ring or the cylinder with bottle attachment can easily be replaced when required.

All components must be thoroughly cleaned and disinfected before each use. Please always follow the instructions for use provided by the respective manufacturers of the cleaning and disinfecting agents used.

1. TO FIT THE NEEDLE

Attention !

- Place the needle onto the Luer-Lock connector and secure it by turning it clockwise.
- Remove the protective cap from the needle immediately before use.
- To remove the needle, first replace the protective cap and then unscrew the needle by turning it counterclockwise.

2. CONNECT THE BOTTLE


1. Before inserting a new bottle, pierce the rubber stopper with a sterile needle.
2. Insert the bottle into the plastic basket from above and press it down evenly until it sits firmly and securely in the basket.
3. Insert the loaded basket into the designated holder on the syringe. Then guide the easyMount basket over the locking mechanism and turn it clockwise until it audibly or noticeably clicks into place. This secures the bottle correctly. The suction tip must be inserted through the centre of the opening into the bottle and sit completely within the bottle attachment.


3. EXCHANGE THE BOTTLE

- Unlock the easyMount basket by turning it counterclockwise and pull it out along with the bottle. Push the bottle out of the basket, replace it, and insert a new bottle.

4. DOSE ADJUSTMENT

- Turn the dosing sleeve clockwise or counterclockwise until the desired dosage is displayed on the cylinder scale. The center of the O-ring serves as the reference for adjustment.
- On certain models, the dosage amount is preset and cannot be changed.

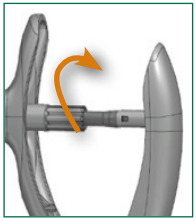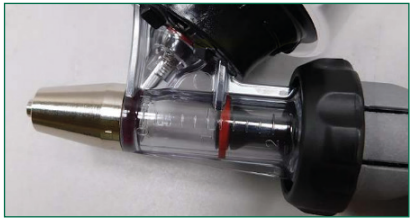Dosing edge

5. BASKET POSITION ADJUSTMENT

1. Loosen the clip nut by about 1.5-2 turns.
2. Gently pull the cylinder forward.
3. Rotate the basket to the desired position to align the bottle correctly.
4. Continue turning the basket clockwise until the bottle is securely locked in place. Please note that the easyMount basket can only be rotated clockwise. Ensure that the cylinder is fully locked in place and can no longer be moved.
5. Then tighten the clip nut securely again.

Different positioning options for the basket

6. FILLING OF THE EM EasyTwist™

1. Remove the protection cover of the needle.
2. Hold the syringe at a slight angle. The suction tip must always remain below the liquid level to prevent air from being drawn in.
3. Pull the plunger back and let it slide forward again to draw liquid into the barrel.
4. If air bubbles are visible, repeat the process until there is no more air in the barrel.
5. When filling, the syringe should be set to at least 0.5 ml. Before use, check that there are no more air bubbles in the cylinder. If necessary, hold the cannula upright to allow trapped air to escape. Afterward, any desired dosage can be set.

INSTRUCTION AFTER USE OF THE EM EasyTwist™

After use, remove all residues with suitable cleaning agents and rinse the EM EasyTwist™ thoroughly with warm water.

Attention !

If the syringe draws air after cleaning, the valve in the front metal attachment may have become stuck. In this case, gently tap the front part of the syringe against a solid object, such as the edge of a table, and then rinse it again with warm water.

DISINFECTION

It is recommended to carefully disinfect the EM EasyTwist™ and the cannulas used before each use.

1. First, clean the syringe with hot water. A bottle of warm water can be used for this.
2. Then place the entire EM EasyTwist™, including the cannulas, in a container of water and boil for about 20 minutes. For a particularly thorough cleaning, it is recommended to separate the cylinder from the handle and boil the syringe in disassembled form.
3. After boiling, remove the syringe, dry it with a clean cloth, and then flush it empty until no residual moisture remains.

TO LUBRICATE

Remove the cylinder by turning it. Then lightly lubricate the piston and cylinder with 100% pure, cold-pressed olive oil or a comparable edible oil. Please follow the oil manufacturer’s storage and shelf-life instructions. Do not use mineral oils, silicone oils, WD-40, glycerin, or alcohol. Petroleumbased oils can damage the O-rings or reduce the product’s lifespan.

© copyright Healthcare FOR Animals GmbH – Subject to technical and design changes without notice

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© copyright Healthcare FOR Animals GmbH – Änderungen in Technik und Design vorbehalten
© copyright Healthcare FOR Animals GmbH – Reservado el derecho a realizar modificaciones técnicas y de diseño
© copyright Healthcare FOR Animals GmbH – Sous réserve de modifications techniques et de conception