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Instruction manual

Manual de instrucciones

Mode d'emploi

# MM Repeater™

**by Healthcare FOR Animals**

**for veterinary use only**

**Healthcare FOR Animals GmbH**
Max-Planck-Straße 14
78549 Spaichingen
Deutschland
info@healthcareforanimals.de

**MM Repeater™**
Von HC4A

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
HEALTHCARE FOR ANIMALS

➢ Dosing syringe for Veterinary use

A) Metal needle hub, Luer Lock (LL)
B) Piston with O-Ring
C) High-strength plastic barrel
D) Dose adjustment wheel
E) Locking button
F) Ergonomic handle
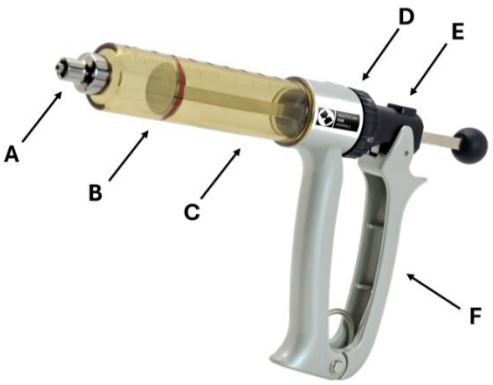
## Instruction for use

The MM Repeater™ repeating syringe has been designed for everyday use under demanding field conditions. To ensure reliable performance and proper hygiene, the device should be thoroughly cleaned and disinfected after each use. Always follow the instructions provided by the manufacturer of the cleaning and disinfecting agents used.

### 1. TO FIT THE NEEDLE

**Attention**: For single-use needles, remove the protective cap only immediately before use.

- Place the needle onto the Luer Lock connection.
- Secure it by turning clockwise until firmly attached.
- To remove the needle, turn it counterclockwise.
- After use, replace the protective cap on the needle.

### 2. FILLING THE SYRINGE

- Hold the syringe by the fixed front section of the handle.
- Press the release button to allow free movement of the plunger rod.
- Pull back the plunger rod and draw up the desired liquid.
- Actuate the movable handle until the plunger rod re-engages.
- The syringe is now ready for use.

### 3. DOSE ADJUSTMENT

- Depending on the model, the required dose can be selected using the adjustment wheel:
- In 1-ml steps for the 50-ml model: 1 – 2 – 3 – 4 – 5 ml
- In 0.5-ml steps for the 25-ml version: 0.5 – 1 – 1.5 – 2 – 2.5 ml

### INSTRUCTION AFTER USE OF THE MM Repeater™

- Remove all product residues after use.
- Clean the syringe using a suitable cleaning agent.
- Thoroughly rinse with warm water.

## DISINFECTION

To ensure hygienic operation, it is recommended to disinfect both the syringe and the needles before and after each use.

1. First clean the syringe with hot water.
2. Boil the syringe, including the attached needle, for up to 20 minutes.
3. Alternatively, the barrel with the metal needle fitting may be autoclaved at 134°C (273°F).
4. After disinfection, remove all components, dry them thoroughly, and ensure any residual moisture is removed from the syringe.

## TO LUBRICATE

Remove the cylinder by turning it counterclockwise. Then lightly lubricate the piston and O-Ring with 100% pure, cold-pressed olive oil or a comparable edible oil.

Please follow the oil manufacturer's storage and shelf-life instructions.

Do not use mineral oils, silicone oils, WD-40, glycerin, or alcohol. Petroleum-based oils can damage the O-rings or reduce the product's lifespan.

© Copyright Healthcare FOR Animals GmbH – Subject to technical and design modifications without notice.

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© copyright Healthcare FOR Animals GmbH – Änderungen in Technik und Design vorbehalten.
© Copyright Healthcare FOR Animals GmbH – Reservado el derecho a realizar modificaciones técnicas y de diseño.
© Copyright Healthcare FOR Animals GmbH – Sous réserve de modifications techniques et de conception.